

# Overbed Table

## PRODUCT MANUAL



**Scan Here**  
With Your Phone To  
**Get Started!**



OVERBED TABLE NON-TILT

 **journey**  
HEALTH & LIFESTYLE®

[www.shopjourney.com](http://www.shopjourney.com)

# Overbed Table

- 1. Safety Warnings ..... 3
- 2. Assembly Instruction ..... 4-5
- 3. How to Use ..... 6
- 4. Care and Maintenance..... 6
- 5. Warranty ..... 7

Revision Date: August 2026

### WARNING

Read all instructions carefully before assembling to avoid possible injury. Journey Health & Lifestyle shall have no liability for any injury resulting from improper assembly or misuse of this product.

- DO NOT overtighten the bolts, screws or use any power tools when assembling this product.
- DO NOT use the product if it appears to be damaged or unstable.
- DO NOT modify the product in any way.
- DO NOT overload the product as that can cause it to break and/or serious injury.
- DO NOT use if transfer bench is damaged. Check frequently for cracks, damage or looseness.
- Maximum weight capacity is 50 lbs.




**Capacity not to exceed 22 lbs**

### STEP 1

Lay tabletop upside down on a flat surface with the pre-drilled holes facing up. Align the mounting bracket to the pre-drilled holes of the tabletop and hand tighten all eight screws.

Use the Allen wrench to securely tighten the screws (see Figure 1).

*Figure 1*

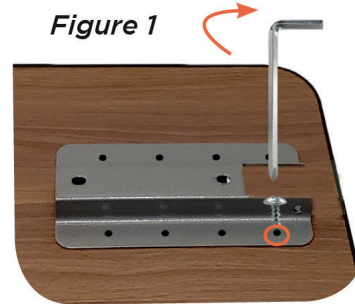

### STEP 2

Stand the column assembly right side up. Be sure the rod is not protruding out of the top. Insert the lever with the rounded side facing up into the column assembly. Ensure that the flat side of the lever rests on the rod (see Figure 2).


*Figure 2*



### STEP 3


While holding the lever securely in place on top of the rod, flip the column assembly upside down and insert the column assembly into the tabletop mounting bracket (see Figure 3).

*Figure 3*

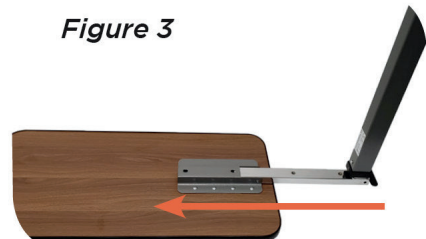

### STEP 4

Align the horizontal hole of the column to the tabletop mounting bracket, and align the holes on the column assembly to the hole through the lever. Insert the lever screw and fully tighten the nut. 10mm wrench size (see Figure 4).

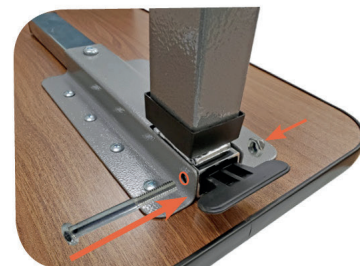


*Figure 4*

## 2. ASSEMBLY INSTRUCTIONS

# Overbed Table

### PARTS INCLUDED

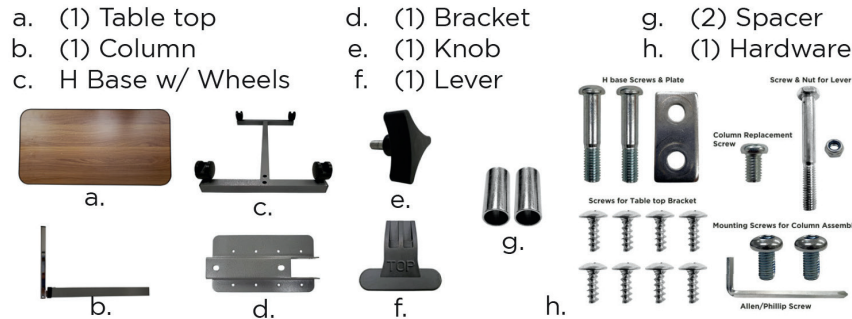

### STEP 5

Secure the column assembly to the mounting bracket with two mounting screws through the vertical holes in the mounting bracket and into the threaded holes in the column assembly (see Figure 5).

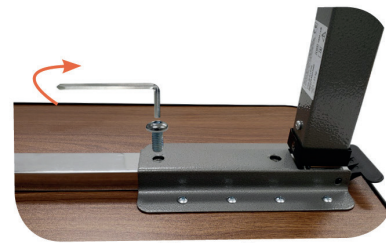


Figure 5

### STEP 6

Locate the two holes on the H base. Insert both bolt sleeves into the bottom of the H base (see Figure 8). Align the H base to the column with wheels facing up while tabletop is upside down. Place the metal plate to the H base holes and tighten both bolts (see Figure 9).

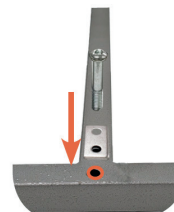


Figure 8



Figure 9

### STEP 7

Turn the table over. Replace the top black screw that's located on the column with the chrome screw that came with the hardware.

**Note: Do not take off black screw until steps 1 - 6 are completed.**



Figure 10



Figure 11

### STEP 8

Insert the knob onto the hole that is located on the side of the column (see Figure 12).

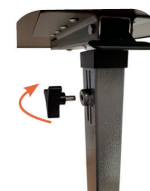


Figure 12

## How to Use

Adjust the height of the table by removing or loosening the knob.

1. To raise - Pull up from the bottom of the tabletop until you reach your desired height.
2. To lower - Pull up on the lever and push down on the tabletop. Release the lever at your desired height.

## Care and Maintenance

Check to make sure all hardware is securely tightened before using the table. Use a mild detergent and warm water to clean. Dry with a soft cloth to avoid moisture buildup.



### **Overbed Table Non-Tilt 12-Month Limited Warranty**

Journey Health & Lifestyle warrants the frame and mechanical components to be free from defects in materials and workmanship for a period of twelve (12) months from the original purchase date. The originally supplied table is all warranted to be free from defects in materials and workmanship for a period of six (6) months from the original purchase date. If within such warranty period the components or parts referenced in this paragraph shall be proven to Journey Health & Lifestyle's satisfaction to be defective, the product shall, at Journey Health & Lifestyle's option, be repaired or replaced with new or refurbished product components. In some cases, an installation kit with instructions may be provided to the owner to facilitate installation of a replacement product component. This warranty does not include any labor charges incurred in product component installation. The warranty period of the repaired or replaced product shall terminate with the termination of the warranty period for the original product. Journey Health & Lifestyle's sole obligation and your exclusive remedy under this warranty shall be limited to such repaired or replaced product.

Contact Journey Health & Lifestyle before sending any product for warranty repair or replacement. Do not send any product without our prior written consent.

### **Limitations and Exclusions:**

This warranty does not cover damage, deterioration or malfunction resulting from accident, negligence, misuse, abuse, occupancy by more than one user and cargo weight above that specified in this Product Manual; improper storage; improper operation or failure to follow instruction according to this Product Manual; damage in shipment of the product (claims must be presented to the carrier); repair, attempted repair or modifications made by anyone other than Journey Health & Lifestyle; any unit on which the serial number has been defaced, modified or removed; normal wear items; any periodic maintenance; deterioration due to perspiration, corrosive atmosphere, or ground covering such as pine needles or crushed limestone dust, or other external causes such as punctures, extremes in temperature or humidity; damages related to lighting, flood or circumstances beyond Journey Health & Lifestyle's control; or product subjected to commercial or institution use.

This warranty applies in North America only and is extended only to the original purchaser of the Overbed Table Non-Tilt product when purchased directly from Journey Health & Lifestyle (USA) or an authorized reseller. The warranty is nontransferable beyond the initial purchaser. This warranty gives you specific legal rights and you may also have other legal rights which vary from place to place.

The foregoing warranty is exclusive and in lieu of all other express warranties, implied warranties, including but not limited to the implied warranties of merchantability and fitness for a particular purpose, and shall not extend beyond the duration of the express warranty provided herein. Journey Health & Lifestyle shall not be liable for any consequential or incidental damages whatsoever and disclaims any liability beyond the initial purchase price of the Overbed Table Non-Tilt product.

Copyright © 2026 Journey Health & Lifestyle

**IMPORTANT WARNING**  
**FOLLOW ALL DIRECTIONS CAREFULLY BEFORE USE.**  
**FAILURE TO DO SO MAY RESULT IN SERIOUS INJURY.**



# Overbed Table

If you have any questions, call us at our toll free number:

**1-800-958-8324**