

Overbed Table

PRODUCT MANUAL



OVERBED TABLE



www.shopjourney.com

HOW TO ASSEMBLE YOUR



OVERBED TABLE

Capacity not to exceed 22 lbs



STEP ONE

Remove the product from its packaging.


STEP TWO

Pull the positioning pin located on the bottom frame towards the depicted figure. Rotate the product vertically until the pin aligns and secures it in place.



STEP THREE

Pull the positioning pins on the table's positioning iron plate in the direction indicated in the diagram. Adjust the table's angle as needed so that the pins align and secure the table.

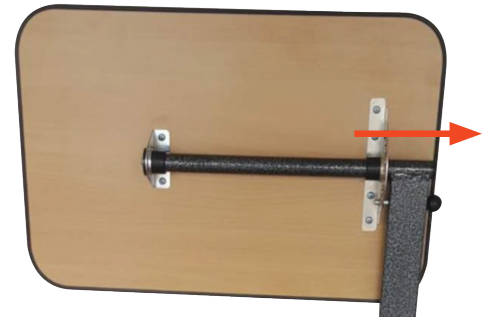

When adjusting the table height, simply pull it up or down to the desired level and secure it using the knob on the tube.

If the table height is too high, lower it by pulling out the positioning pin on the pipe fitting. Adjust to the desired level and release the positioning pin to hold it in place.



Read all instructions carefully before assembling to avoid possible injury. Journey Health & Lifestyle shall have no liability for any injury resulting from improper assembly or misuse of this product.

- DO NOT overtighten the bolts, screws or use any powertools when assembling this product.
- DO NOT use the product if it appears to be damaged or unstable.
- DO NOT modify the product in any way.
- DO NOT overload the product as that can cause it to break and/or serious injury.
- DO NOT use if transfer bench is damaged. Check frequently for cracks, damage or looseness.
- Maximum weight capacity is 22 lbs.

4. HOW TO USE

Adjust the height of the table by removing or loosening the knob.

1. **To Raise** — Pull up from the bottom of the tabletop until you reach your desired height.
2. **To Lower** — Pull up on the lever and push down on the tabletop. Release the lever at your desired height.

5. CARE AND MAINTENANCE

Check to make sure all hardware is securely tightened before using the table. Use a mild detergent and warm water to clean. Dry with a soft cloth to avoid moisture buildup.

Overbed Table 12-Month Limited Warranty

Journey Health & Lifestyle warrants the frame and mechanical components to be free from defects in materials and workmanship for a period of twelve (12) months from the original purchase date. The originally supplied table is all warranted to be free from defects in materials and workmanship for a period of six (6) months from the original purchase date. If within such warranty period the components or parts referenced in this paragraph shall be proven to Journey Health & Lifestyle's satisfaction to be defective, the product shall, at Journey Health & Lifestyle's option, be repaired or replaced with new or refurbished product components. In some cases, an installation kit with instructions may be provided to the owner to facilitate installation of a replacement product component. This warranty does not include any labor charges incurred in product component installation. The warranty period of the repaired or replaced product shall terminate with the termination of the warranty period for the original product. Journey Health & Lifestyle's sole obligation and your exclusive remedy under this warranty shall be limited to such repaired or replaced product.

Contact Journey Health & Lifestyle before sending any product for warranty repair or replacement. Do not send any product without our prior written consent.

Limitations and Exclusions:

This warranty does not cover damage, deterioration or malfunction resulting from accident, negligence, misuse, abuse, occupancy by more than one user and cargo weight above that specified in this Product Manual; improper storage; improper operation or failure to follow instruction according to this Product Manual; damage in shipment of the product (claims must be presented to the carrier); repair, attempted repair or modifications made by anyone other than Journey Health & Lifestyle; any unit on which the serial number has been defaced, modified or removed; normal wear items; any periodic maintenance; deterioration due to perspiration, corrosive atmosphere, or ground covering such as pine needles or crushed limestone dust, or other external causes such as punctures, extremes in temperature or humidity; damages related to lighting, flood or circumstances beyond Journey Health & Lifestyle's control; or product subjected to commercial or institution use.

This warranty applies in North America only and is extended only to the original purchaser of the Overbed Table product when purchased directly from Journey Health & Lifestyle (USA) or an authorized reseller. The warranty is nontransferable beyond the initial purchaser. This warranty gives you specific legal rights and you may also have other legal rights which vary from place to place.

The foregoing warranty is exclusive and in lieu of all other express warranties, implied warranties, including but not limited to the implied warranties of merchantability and fitness for a particular purpose, and shall not extend beyond the duration of the express warranty provided herein. Journey Health & Lifestyle shall not be liable for any consequential or incidental damages whatsoever and disclaims any liability beyond the initial purchase price of the Overbed Table product.

Copyright © 2024 Journey Health & Lifestyle

IMPORTANT WARNING
FOLLOW ALL DIRECTIONS CAREFULLY BEFORE USE.
FAILURE TO DO SO MAY RESULT IN SERIOUS INJURY.

- Manufacturer: Combo
- Product Number: 08982
- Tabletop Adjustment: 700-930mm 35mm
- Tray Length: 600mm
- Tray Width 400mm
- Base Length: 58635mm
- Weight: 12.0KGS
- Weight Limit:
- Tabletop Weight Capacity 10KG
- Warranty: Lifetime Limited 1 Year

Overbed Table



Overbed Table

If you have any questions, call us at our toll free number:

1-800-958-8324