

Overbed Table

PRODUCT MANUAL



Scan Here
With Your Phone To
Get Started!



OVERBED TABLE

 **journey**
HEALTH & LIFESTYLE®

www.shopjourney.com

- 1. Introduction 3
- 2. Assembly Instruction 4-5
- 3. How to Use 6
- 4. Care and Maintenance..... 6
- 5. Warranty 7
- 6. Technical Specifications..... 8

Welcome to Journey Health & Lifestyle

Thank you for purchasing your foldable Overbed Table. It's designed to bring everyday convenience, comfort, and independence right to the bedside. Whether you're recovering, working from bed, enjoying meals, or reading comfortably, this table adapts to your needs.

With smooth height adjustment, a stable rolling base, and a wide, easy-clean tabletop, it creates a dependable surface for daily activities without strain or awkward positioning. The low-profile frame slides effortlessly under most beds, recliners, and chairs, allowing the table to position exactly where you need it: close, secure, and accessible.

Safety Warnings

- Read all instructions carefully before assembly and use to avoid possible injury.
- Do not overtighten bolts or screws, and do not use power tools during assembly.
- Do not modify the product in any way.
- Do not overload the table; exceeding the 22 lb capacity may cause damage or injury.
- Do not lean, sit, or stand on the tabletop.
- Place heavier items toward the center of the tabletop, not at the far edges, to reduce tipping risk.
- Do not use the table as a rail or barrier, and do not rely on it to prevent falls from the bed. People who try to climb over or around it are at increased risk of injury.
- People who are confused, very frail, or have dementia may be at higher risk of getting caught between the table and the bed; check positioning carefully and supervise when needed.
- Check that all knobs, pins, and height-adjustment parts are fully engaged and tightened before placing items on the table or moving it.
- When adjusting height or angle, keep fingers clear of moving joints and hinges to avoid pinching.
- Move the table slowly and carefully, especially over uneven flooring or thresholds, to prevent tipping.
- Do not use the table if any part appears bent, cracked, or unstable; contact the supplier before continued use.

Important Capacity Limit

The maximum load for this overbed table is 22 lbs. Do not place items heavier than this on the tabletop.

Assembly Instructions

Step 1: Unpack the Table

Remove the overbed table and all packaging materials from the box. Place the table on a flat, stable surface for assembly.


Step 2: Set the Base in Position

Locate the positioning pin on the bottom frame. Pull the pin toward the direction shown in the diagram, then rotate the frame vertically until the pin aligns with its hole and locks the frame securely in place.



Step 3: Set the Tabletop Angle

On the table's positioning iron plate, pull the positioning pins in the direction indicated in the diagram. Adjust the tabletop angle to the position you prefer, then release the pins so they align with the holes and secure the table at that angle.

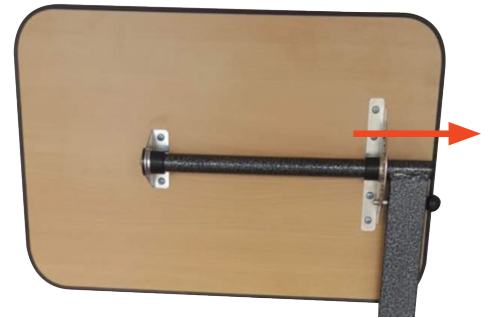

Height Adjustment

- To raise or lower the table, slide the tabletop up or down to the desired height, then tighten the knob on the support tube to hold it in place.
- If the table is set too high, pull out the positioning pin on the pipe fitting, lower the table to a comfortable height, and release the pin to lock it.



How to Use

- To raise the table: Loosen or remove the height adjustment knob if needed, then pull up on the tabletop from underneath until you reach your desired height and retighten the knob.
- To lower the table: Pull up on the height adjustment lever and gently push down on the tabletop, then release the lever at the desired height.
- Before each use, check that all hardware is securely tightened and the table is stable.

Care and Maintenance

Clean the table with a mild detergent and warm water, then dry with a soft cloth to prevent moisture buildup. Periodically check that all fasteners are tight and that the table moves and locks smoothly at different heights.

Overbed Table 12-Month Limited Warranty

Journey Health & Lifestyle warrants the frame and mechanical components to be free from defects in materials and workmanship for a period of twelve (12) months from the original purchase date. The originally supplied table is all warranted to be free from defects in materials and workmanship for a period of six (6) months from the original purchase date. If within such warranty period the components or parts referenced in this paragraph shall be proven to Journey Health & Lifestyle's satisfaction to be defective, the product shall, at Journey Health & Lifestyle's option, be repaired or replaced with new or refurbished product components. In some cases, an installation kit with instructions may be provided to the owner to facilitate installation of a replacement product component. This warranty does not include any labor charges incurred in product component installation. The warranty period of the repaired or replaced product shall terminate with the termination of the warranty period for the original product. Journey Health & Lifestyle's sole obligation and your exclusive remedy under this warranty shall be limited to such repaired or replaced product.

Contact Journey Health & Lifestyle before sending any product for warranty repair or replacement. Do not send any product without our prior written consent.

Limitations and Exclusions:

This warranty does not cover damage, deterioration or malfunction resulting from accident, negligence, misuse, abuse, occupancy by more than one user and cargo weight above that specified in this Product Manual; improper storage; improper operation or failure to follow instruction according to this Product Manual; damage in shipment of the product (claims must be presented to the carrier); repair, attempted repair or modifications made by anyone other than Journey Health & Lifestyle; any unit on which the serial number has been defaced, modified or removed; normal wear items; any periodic maintenance; deterioration due to perspiration, corrosive atmosphere, or ground covering such as pine needles or crushed limestone dust, or other external causes such as punctures, extremes in temperature or humidity; damages related to lighting, flood or circumstances beyond Journey Health & Lifestyle's control; or product subjected to commercial or institution use.

This warranty applies in North America only and is extended only to the original purchaser of the Overbed Table product when purchased directly from Journey Health & Lifestyle (USA) or an authorized reseller. The warranty is nontransferable beyond the initial purchaser. This warranty gives you specific legal rights and you may also have other legal rights which vary from place to place.

The foregoing warranty is exclusive and in lieu of all other express warranties, implied warranties, including but not limited to the implied warranties of merchantability and fitness for a particular purpose, and shall not extend beyond the duration of the express warranty provided herein. Journey Health & Lifestyle shall not be liable for any consequential or incidental damages whatsoever and disclaims any liability beyond the initial purchase price of the Overbed Table product.

Copyright © 2026 Journey Health & Lifestyle

IMPORTANT WARNING
FOLLOW ALL DIRECTIONS CAREFULLY BEFORE USE.
FAILURE TO DO SO MAY RESULT IN SERIOUS INJURY.

- Manufacturer: Combo
- Product Number: 08982
- Tabletop Adjustment: 700-930mm 35mm
- Tray Length: 600mm
- Tray Width 400mm
- Base Length: 58635mm
- Weight: 12.0KGS
- Weight Limit:
- Tabletop Weight Capacity 10KG
- Warranty: Lifetime Limited 1 Year



Overbed Table

If you have any questions, call us at our toll free number:

1-800-958-8324